VibroSculpt™

INSTRUCTIONS MANUAL

MANUAL DE INSTRUCCIONES

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MANUEL D'INSTRUCTIONS

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EN Features:

- Portable and easy to use.
- Interchangeable pads.
- Variable speed for greater comfort.
- Ergonomic design.

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EN ERGONOMIC MASSAGER
The ergonomically designed
Vibro Sculpt™ massage device is
the new portable massager
designed to be used in the
comfort of your home.

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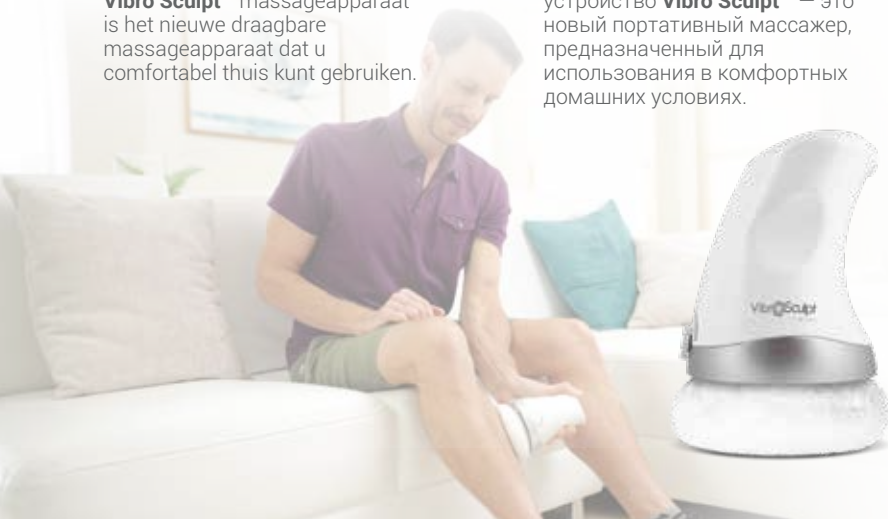
EN VIBROBUFF TECHNOLOGY
Vibro Sculpt™ is very easy to use; by following the instructions, it will soon become part of your daily routine.

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EN Benefits:

- Gentle and pleasant massage.
- Ideal for daily routine.
- Offers an experience of self-care and personal well-being.
- Enjoy a relaxing moment anywhere.

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1

EN The **Vibro Sculpt™** massager is designed to be used **ONLY** with the exclusive Beauty Pro pads.

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2

EN Make sure the switch is in the **OFF** position.

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3

EN Plug the adapter into the wall outlet and ensure the voltage matches the mains voltage.

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4

EN Choose the massage pad for the desired area.

BLUE: for a gentle massage

WHITE: for a relaxing massage

ORANGE: for use in combination with body gels and creams

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EN HOW TO MASSAGE THE ABDOMEN AND WAIST

- Lie down in bed or stand in a comfortable position. Keep your back straight and your knees slightly bent. Starting from your navel, apply the device at a low level over your entire abdomen, then move it to your sides. Leave it on for a couple of minutes on each side of your waist.
- Or you can start from the middle of your abdomen and move outward in straight lines to the sides. Repeat this motion several times, choosing a different starting point each time.
- This massage should be brief, never more than 5-15 minutes.

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EN HOW TO MASSAGE YOUR LEGS

- To massage your legs, use it in upward motions. Sit or lie down in a semi-reclining position and let the device work for a couple of minutes on the inner and outer sides of your knees.
- Continue on both sides of your thighs until you reach your groin and hips.
- Leave it on for 5-15 minutes, always moving slowly upward.
- Hold it firmly, but without excessive pressure. Move slowly and decisively, avoiding quick or imprecise movements.

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EN HOW TO MASSAGE YOUR GLUTES AND HIPS

- Use this device from the base of your buttocks to your waist, passing through your hips.
- The best position is to lie on your back, with your knees bent to slightly raise your hips and leave room for the massage.
- You can also sit or lie down semi-reclined, stretch your legs, and turn your body to one side, resting on your hips, to massage one buttock. When you're finished, switch sides and massage the other.

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EN HOW TO MASSAGE YOUR BACK

- Sit in bed, on the couch, or stand in a comfortable position, keeping your back straight. Place the device on the center of your back and slide it sideways on the low setting, leaving it for a couple of minutes on each side of your back.
- Another option is to start from the center of your back and move it sideways in straight lines. Repeat this motion several times, choosing a different starting point each time.
- This massage should be brief, lasting 5-15 minutes.

DO NOT USE ON THE NECK AREA AND KEEP IT AWAY FROM LONG HAIR, JEWELRY, AND LOOSE CLOTHING.

DO NOT USE ON THE HEAD, FACE, NECK, CERVICAL SPINE, OR BREASTS.

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EN SAFETY INFORMATION

Keep the unit and its packaging materials out of the reach of children.

- Only work with the supplied power supply. If there is obvious damage to the device or its accessories, do not use the device.
- Keep the unit, power supply, and cable away from water, moisture, or other liquids.
- Never use it in a bathroom, sauna, etc.
- Never attempt to retrieve the device if it falls into water. Immediately disconnect the plug from the electrical outlet.
- Repairs to electrical appliances should only be performed by the manufacturer or qualified personnel. Improper repairs can pose a considerable danger to the user.
- Turn off the device and disconnect the power supply after each use and before each cleaning.

Follow local regulations regarding disposal of materials. If you have any questions, please contact the relevant authorities responsible for waste disposal in your area.

PROPER USE

This device is intended only as a massager. The device is for home use only and is NOT intended for medical or commercial purposes. It cannot replace medical treatment. DO NOT use the massage device if any of the following warnings apply to you. Consult your doctor if you are unsure whether the massage device is right for you. DO NOT use the device directly on varicose veins or after surgery. We always recommend consulting your doctor before using this device.

RETURN POLICY

This product is guaranteed by the manufacturer. We will only accept the return of an unpaid, unused Vibro Sculpt™ within 30 days of purchase. Once the seal is broken, the Vibro Sculpt™ device is non-returnable. The box must be in its original condition and must NOT be defaced in any way.

WARNING

Consult your doctor before using this product. Read the manual carefully. Use this product only for its intended use as described in the manual. This product is NOT intended as a substitute for professional medical advice, diagnosis, or treatment. Use this product at your own risk. The manufacturer and distributor assume no responsibility for the use of this product. Close supervision is necessary when this product is used near children or on or by disabled persons. Do not use this product when in pain or if you are pregnant. Do not massage any area of the body that is swollen, inflamed, or covered with a rash.

ALL WARRANTIES, WHETHER EXPRESS OR IMPLIED, OF MERCHANTABILITY, OR FITNESS FOR A PARTICULAR

PURPOSE, ARE HEREBY DISCLAIMED. THE MANUFACTURER OR DISTRIBUTOR SHALL NOT BE LIABLE FOR ANY DAMAGES OF ANY KIND ARISING FROM THE USE OF THIS PRODUCT, INCLUDING, BUT NOT LIMITED TO, DIRECT, INDIRECT, INCIDENTAL, PUNITIVE, AND CONSEQUENTIAL DAMAGES. THE MANUFACTURER'S OR DISTRIBUTOR'S MAXIMUM LIABILITY TO YOU, AND YOUR SOLE AND EXCLUSIVE REMEDY FOR THE USE OF THIS PRODUCT, IS THE AMOUNT PAID BY YOU TO THE MANUFACTURER OR DISTRIBUTOR.

ATTENTION - DISPOSAL OF MATERIALS

A crossed-out wheeled bin symbol means you should inform yourself and follow local regulations regarding the disposal of this type of product.

Do not dispose of this product as you would other household waste.

Dispose of your device in accordance with applicable local regulations.

Electrical and electronic devices contain hazardous substances that may have harmful effects on the environment and/or human health and must be recycled properly.

PRECAUTIONS:

- Consult your physician before using the massager if you suffer from any medical condition, are undergoing any treatment, or are pregnant.
- This is NOT a medical device, and no medical claims are implied or guaranteed by the use of this product.
- Confirm that the unit is turned off before connecting the power plug.
- It is recommended to start the massager at a low speed/intensity first and then gradually increase the speed/intensity as needed.
- DO NOT massage a particular area for more than 15 minutes.
- DO NOT leave the product where children can play with it.
- DO NOT use on varicose veins, moles, or damaged skin tissue.
- DO NOT use on any visible or recent scars.
- Avoid going above the navel. If the abdomen is hard and hot, do not use the device and consult your doctor.
- You may notice some temporary redness on the skin you have treated normally.
- If the product appears to be defective, DO NOT attempt to open or repair it. Repairs should only be performed by the manufacturer or authorized distributor.
- WARNING: NOT FOR USE BY PERSONS UNDER 18 YEARS OF AGE. DO NOT USE IF PREGNANT OR BREASTFEEDING. As with any massage product,

individuals with special medical conditions are advised to consult their physician before use.

- DO NOT USE ON THE NECK AREA AND KEEP IT CLEAN OF LONG HAIR, JEWELRY, AND LOOSE CLOTHING.
- DO NOT USE ON THE HEAD, FACE, NECK, CERVICAL SPINE, OR BREASTS.

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